


ENV 48

ENVIRONMENTAL

Analysis of Herbicides on BPX5

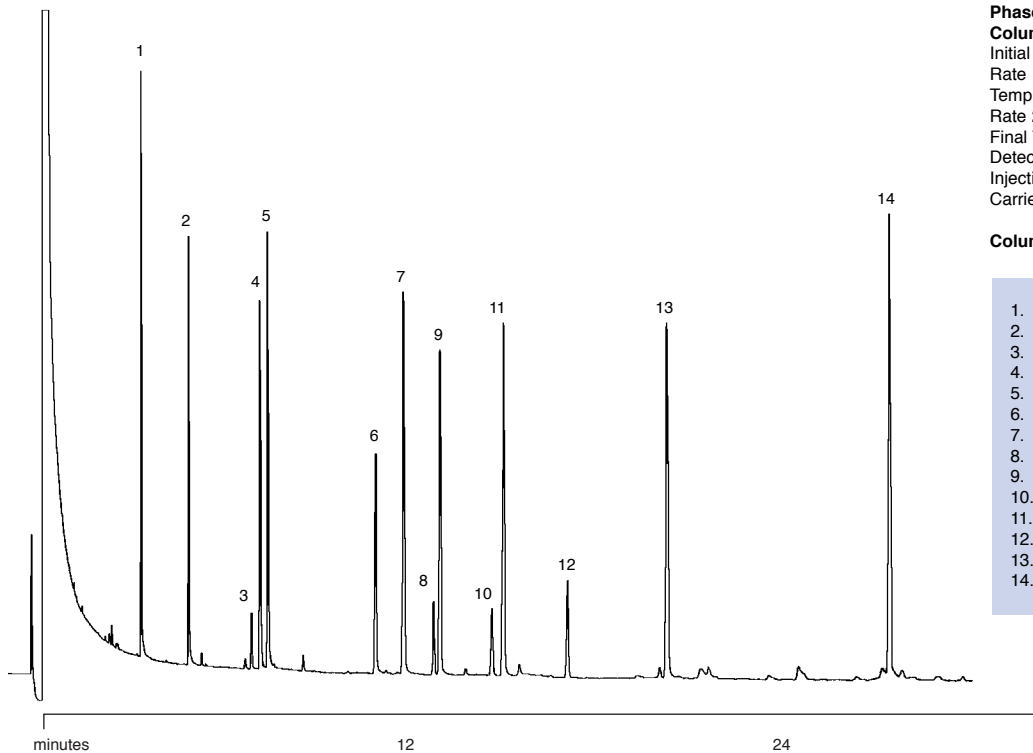

Phase: BPX5, 0.25 μ m
Column: 30m x 0.25mm ID
Initial Temp.: 90°C, 1min
Rate 1: 30°C/min
Temp.: 180°C
Rate 2: 5°C/min
Final Temp.: 260°C, 10min
Detector: NPD
Injection Mode: Varian SPI
Carrier Gas: He, 10psi

Column Part No. 054101

1. Metamidofos
2. Acephate
3. Diphenylamine
4. Monocrofos
5. Sulfotep
6. Tolclofos-methyl
7. Fenitrothion
8. Triadimefon
9. Trichloronate
10. Triadimenol
11. Bromophos-ethyl
12. Bupirimate
13. Carbophenothion
14. Dialifos

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in